


Lower Profile is better

Mini-Mono® Brackets (Low Profile)



Bracket wings deeper undercut

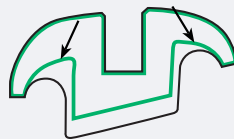

Sandblasted mesh base



Mc Laughlin/Bennett/Trevisi* S₂ System



S-Profil



Low Profile S₂

- **Lower profile** for improved patient comfort
- Mesio-distal width similar to standard bracket for better **rotational control**

 **FORESTADENT**®

Bernhard Förster GmbH · Westliche Karl-Friedrich-Straße 151 · 75172 Pforzheim · Germany
Phone + 49 7231 459-0 · Fax + 49 7231 459-102 · info@forestadent.com · www.forestadent.com


Forestadent USA · 2301 Weldon Parkway · St. Louis, MO 63146 · TOLL FREE 1-800-721-4940
Fax 314-878-7604 · info@forestadentusa.com · Internet: <http://www.forestadentusa.com>

*The FORESTADENT version of these prescriptions are not claimed to be a duplication of any other, nor does Forestadent imply that it is endorsed in any way by the doctor.

Mini-Mono® Brackets (Low Profile)



Lower profile height of the bracket leads to less irritation of the labial mucosa and better aesthetic look.

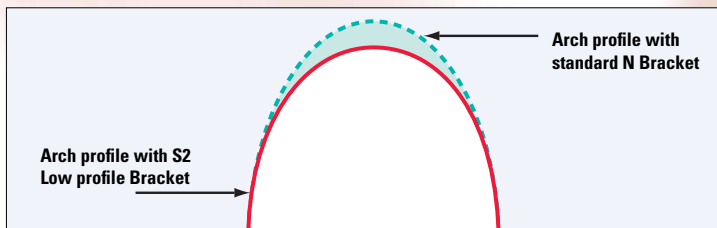

Reduced bracket height changes archwire profile. Standard archwires may be used. Distal cuts will be somewhat larger.




Good rotational control

The mesio-distal width is similar to standard brackets for a better rotational control of the low profile bracket. The occluso-gingival bracket size is smaller and more aesthetically pleasing compared to standard brackets, which also improves patient hygiene.




Bracket wings with deeper undercut for **easy ligating – single and double** eg. with ligature and power chain. Small elastic ligatures have been designed for this bracket (Order no. 659-0004).

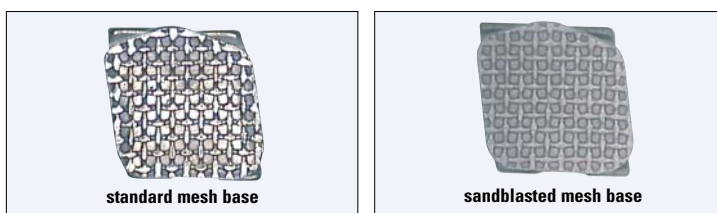

Optimized retention

through sandblasted mesh base.

Slot .022"	Mc Laughlin/Bennett/Trevisi* S ₂ System					CE 0297
Cases	1	5	10	20	50	
Brackets	20	100	200	400	1000	
 Non-Extraction 3**	706-0801	706-0802	706-0803	706-0804	706-0805	
Non-Extraction 3-5**	706-0806	706-0807	706-0808	706-0809	706-0810	
Non-Extraction 3-5**	706-0811	706-0812	706-0813	706-0814	706-0815	
Extraction***	706-0816	706-0817	706-0818	706-0819	706-0820	

**with hook $\frac{3}{3}$ 13 0° torque with hook
***with extraction-brackets $\frac{3}{3}$ 13 0° torque with hook



Bernhard Förster GmbH · Westliche Karl-Friedrich-Straße 151 · 75172 Pforzheim · Germany
Phone + 49 7231 459-0 · Fax + 49 7231 459-102 · info@forestadent.com · www.forestadent.com
Forestadent USA · 2301 Weldon Parkway · St. Louis, MO 63146 · TOLL FREE 1-800-721-4940
Fax 314-878-7604 · info@forestadentusa.com · Internet: http://www.forestadentusa.com

*The FORESTADENT version of these prescriptions are not claimed to be a duplication of any other, nor does Forestadent imply that it is endorsed in any way by the doctor.